
## Advanced Teaching Compound Microscope Hochleistungs-Kursmikroskop

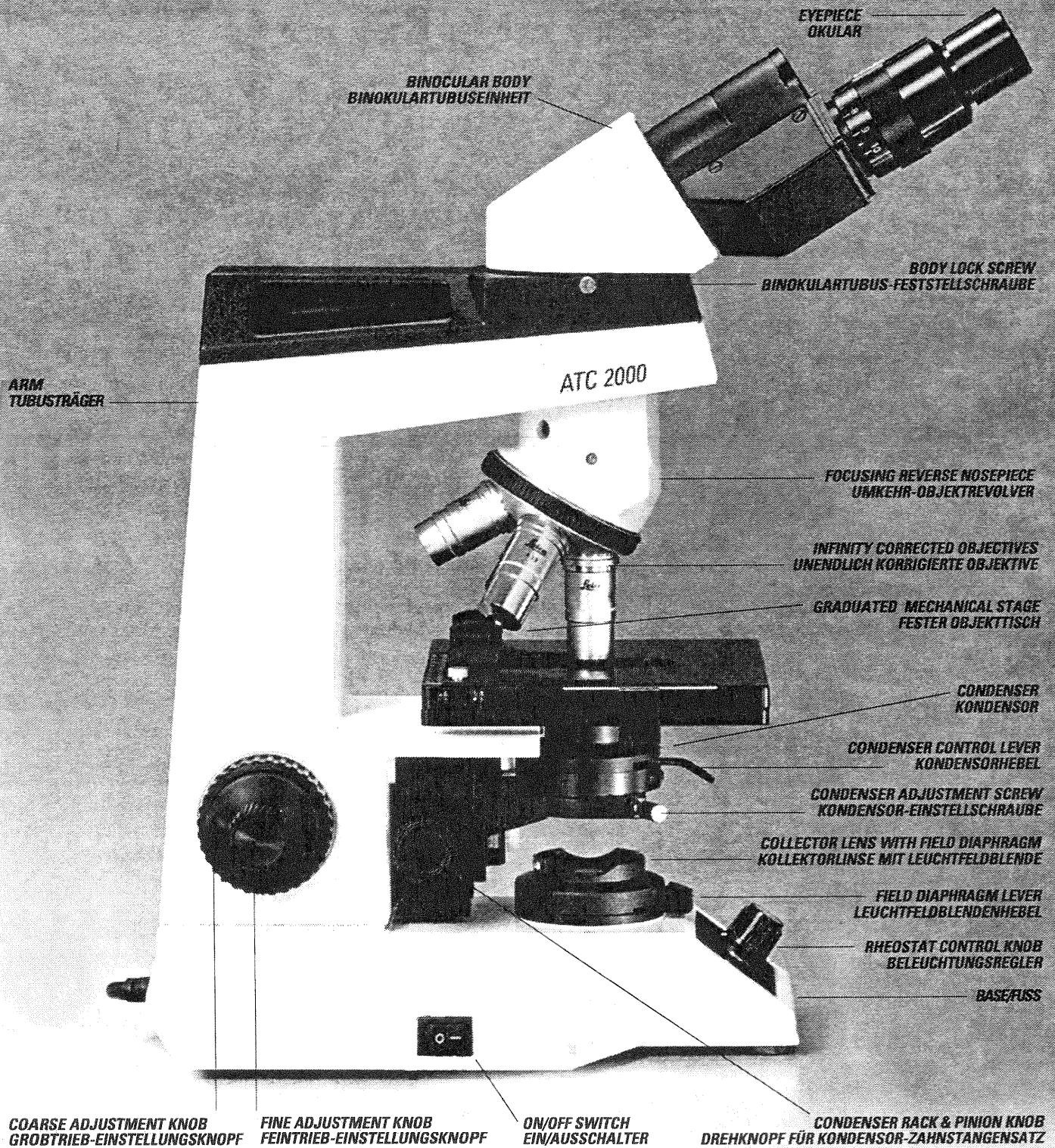**LEICA ATC™ 2000**  
**Advanced Teaching Compound Microscope**  
**Instruction Manual**

**Table of Contents**

<b>Section</b>	<b>Page</b>
1.0 INTRODUCTION .....	2
2.0 UNPACKING .....	2
3.0 PREPARATION FOR USE .....	3
3.1 Infinity Corrected Objectives .....	3
3.2 Viewing Body .....	3
3.3 Substage Condenser .....	3
3.4 Filters .....	3
4.0 OPERATION .....	4
5.0 PHASE CONTRAST MICROSCOPY .....	5
6.0 CARE OF THE MICROSCOPE .....	6
6.1 General .....	6
6.2 Mechanical Maintenance .....	6
6.3 Lamp Replacement .....	7
6.4 Electrical Considerations/Equipment Ratings .....	7
7.0 FUSE REPLACEMENT .....	7

**Explanation of Warning Symbols**

 5X20 Type T200mA Slo-Blo 250V	 <b>CAUTION:</b> Risk of electric shock	 <b>CAUTION:</b> (Refer to accompanying documents)
--	--	---

**1.0 INTRODUCTION**

The LEICA ATC™ 2000 Advanced Teaching Compound Microscope combines exclusive design features and a full range of accessories making it a truly versatile, quality instrument for advanced academic applications.

**2.0 UNPACKING**

Carefully remove the microscope components from the packing carton. Check the components against the following standard equipment list:

1. Stand - includes supporting arm, focusing mechanism, graduated mechanical stage, nosepiece and integral illumination system
2. Viewing Body - monocular, binocular or trinocular, depending upon the component ordered
3. Objectives - as ordered
4. Eyepiece(s) - as ordered
5. Abbe or Phase Contrast Condenser
6. Blue and/or Green Filter
7. Immersion Oil - as ordered
8. Spare Lamp(s) - as ordered
9. Dust Cover

Optional accessory items such as camera adapters and carrying case are not shipped as part of the standard equipment. These items are shipped in separate packages.

Please do not discard the packing material. It can be used for safely storing or transporting the instrument should the need arise.

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**3.1 Infinity Corrected Objectives**

Most microscopes have been constructed with a fixed distance between eyepiece and objective. Leica's infinity corrected system makes it possible to change the objective to eyepiece distance without affecting optical performance. (See Table 1 - Characteristics of Leica Infinity Corrected Objectives and Eyepieces.) At the same time, it simplifies the focusing mechanism; only the nosepiece is focused, not the heavy stage assembly or body tube. Parallel light rays emerging from the objective are picked up by a lens in the microscope body tube and are focused at the correct position in the eyepiece.

Remove the objective lens units from their plastic vials and attach them to the revolving nosepiece holes counterclockwise from the lowest magnification to the highest. Always rotate the nosepiece by screw firmly without forcing it. This precisely locates the body on to be parfocal and parcentered at any selected interpupillary distance setting.

**3.2 Viewing Body**

Attach the viewing body to the stand by backing off the knurled screw, inserting the dovetail into the body support and tightening the knurled screw firmly without forcing it. This precisely locates the body on the optical axis of the microscope, regardless of the selected orientation.

The binocular and trinocular bodies maintain a constant tube length for all interpupillary settings. This means that a change of interpupillary distance does not affect parfocality, magnification or calibrations which depend on magnification. Interpupillary distance is changed by simply spreading or pushing eyetubes. They will hold their position at the appropriate setting. After the body is attached, install eyepieces in body eyepiece tubes.

The left eyepiece tube is focusable to compensate for refractive differences of the eyes. To correctly set the eyetubes, focus on the specimen through the right eyepiece tube only. Use the fine adjustment knob while covering the left eyepiece or closing the left eye. Next, focus the specimen through the left eyepiece by turning the eyetube. Cover the right eyepiece while doing this and be sure to focus with the left eyetube only, without using the focusing knob.

**3.3 Substage Condenser**

Adjust the substage to its lowest position. Push the substage condenser into the fork mount underneath the stage by aligning the positioning screw into the slot. Two knurled centering screws are included to properly align the substage condenser. Each condenser is furnished with an iris diaphragm which controls the effective numerical aperture of the objective. To open and close this diaphragm, simply slide the condenser lever to the right or to the left. To remove the condenser, lower the substage and loosen the screws. Then grasp the condenser and pull forward.

**3.4 Filters**

The blue filter should be placed in the well of the field diaphragm in the base of the instrument. For phase contrast applications replace the blue filter with the green filter. Special care must be used in handling filters for the LEICA ATC™ 2000. Scratches and finger prints will result in poor images and lower light levels.

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#### 4.0 OPERATION

**CAUTION:** A safety precaution, the power cord has been grounded to the microscope base. Never use a 2-prong adapter between the power cord and the power source; it will render the grounding feature ineffective.

- 4.1 Always use your microscope on a hard, stable surface.
- 4.2 Plug the microscope power cord into a suitable grounded electrical outlet. A grounded 3-wire cord is provided.
- 4.3 Turn on the microscope.
- 4.4 Set the rheostat control knob to the lowest setting. The rheostat control allows you to adjust the light level emitted from the base illuminator.
- 4.5 Fully open the field diaphragm by turning the control lever clockwise. Also fully open the aperture diaphragm of the condenser by moving the lever to the extreme right.
- 4.6 Using the substage condenser focusing knob, raise the condenser to the top of its excursion.
- 4.7 Place a specimen slide on the stage.
- 4.8 Rotate the nosepiece to move the 10X objective into working position.
- 4.9 Lower the 10X objective by rotating the coarse adjustment knob to its positive stop. Using the fine adjustment knob, bring the specimen into sharp focus.
- 4.10 Adjust the eyetubes for interpupillary distance and eye difference.
- 4.11 While viewing through the microscope, partially close the field diaphragm using the field diaphragm lever until the iris leaves are imaged within the field of view. Using the condenser focusing knob, adjust the condenser until the leaves are in sharp focus.
- 4.12 Center the image of the field diaphragm by turning the condenser centering screws. This is best accomplished by rotating the two screws simultaneously. After centering, open the field diaphragm until the iris leaves disappear just beyond the field of view.

**CAUTION:** This instrument should only be used as described. Hazards may result if used improperly.

- 4.13 Remove an eyepiece and view the back aperture of the objective. Close the condenser aperture diaphragm and then, to obtain the full resolving power of the microscope, reopen until the diaphragm leaves disappear from view. Replace the eyepiece. The aperture diaphragm can be adjusted to enhance contrast and/or increase the depth of focus.
- 4.14 When changing to higher power objectives, the positions of the aperture and field diaphragms must be reset. As magnification increases, the aperture diaphragms must be opened and the field diaphragm must be opened and closed as required. Also, recenter the field diaphragm and refocus the condenser, if necessary.

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## 5.0 PHASE CONTRAST MICROSCOPE

Many microscopic objects are seen because they have regions of varying density. In normal brightfield illumination, a completely transparent specimen is very difficult to see in any detail; all parts are equally dense. A form of illumination, known as phase contrast, is of value in the study of transparent media.

Phase contrast is a method of illumination in which a portion of the light is treated differently from the rest. This portion interferes with the remainder to produce a visible image of an otherwise invisible, transparent specimen. Standard phase accessories include 10X, 20X, 40X and 100X phase objectives, phase turret condenser, green filter and telescopic eyepiece. Instructions for assembling phase contrast microscopes are the same as for the standard brightfield microscope.

Once the instrument has been assembled, it is ready for the final alignment for phase contrast.

Rack the substage phase condenser to the top of its travel. Swing the 10X objective into the optical path and rotate the phase turret to align the number 10 with the index mark. Focus the microscope on the specimen in the normal manner. Remove either one of the eyepieces and replace it with the phase telescope. Loosen the small locking screw on the telescope and while looking through the telescope raise or lower the inner tube until the phase ring in the objective is in focus. Now adjust the condenser focusing knob to bring the phase annuli into focus. Note the two annuli. The darkened larger annulus is the phase retardation optic mounted in the microscope objective. The brighter and smaller annulus is an image of the illumination annulus located in the phase condenser.

To achieve optimum performance, the bright annulus must be positioned concentric with the dark annulus. If it is not concentric, it must be adjusted. This can be done by adjusting the two knurled centering screws located on the front right and left of the turret condenser.

When the two annuli are accurately aligned with the 10X objective, align the remaining annuli the same way. Once this is done, remove the phase telescope and reinstall the eyepiece into the microscope body. The microscope is now ready to use.

The green filter is supplied as part of the phase accessories and can be used by placing it into the light port in the microscope base. Simply remove the blue filter from the light port in the microscope base and replace it with the green filter.

The objectives in the phase kit will also perform well as normal brightfield objectives. To change from phase to brightfield viewing, simply rotate the condenser turret to the "0" position and adjust the aperture iris and the field diaphragm for optimum contrast.

Catalog Number	Power	Type	Numerical Aperture	Working Distance (mm)	Field of View (mm) 10X Eyepieces No. 191/191A
438	4X	Achromat	0.10	8.83	5.0
439	10X	Achromat	0.25	0.66	2.0
440	40X	Achromat	0.60	0.49	0.5
467	100X(Oil)	Achromat	1.25	0.12	0.2
464	4X	Planachromat	0.10	6.06	5.0
465A	10X	Planachromat	0.25	2.51	2.0
466	40X	Planachromat	0.65	0.55	0.5
463	100X(Oil)	Planachromat	1.25	0.09	0.2
1750	10X	Phase Planachromat	0.25	4.15	2.0
1752	20X	Phase Planachromat	0.50	1.50	1.0
1751	40X	Phase Planachromat	0.65	0.55	0.5
1753	100X(Oil)	Phase Planachromat	1.25	0.09	0.2

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## **6.0 CARE OF THE MICROSCOPE**

### **6.1 General**

Keeping all optical components clean is important to good optical performance. The microscope should always be covered with the plastic dust cover provided with the instrument when it is not in use. If any optical surface becomes badly coated with dust or dirt, it should be cleaned by blowing it off with a syringe or brushing it off with a camel hair brush before attempting to wipe the surface clean.

Optical surfaces should be cleaned with a lint-free cloth, lens tissue or a cotton swab moistened with methanol or commercially available glass cleaner. It is very important to avoid the excessive use of solvent. The cloth lens tissue or cotton swab should be moistened with solvent, but not be wet enough for the solvent to seep around the lens with the resultant danger of loosening cement on interior surfaces.

No part of the microscope is quite so vulnerable to dirt, dust and oil as the front lens of the objective. Whenever lack of contrast, cloudiness or poor definition is encountered, carefully check the condition of the front lens with a magnifier.

The lower magnification objectives have fairly large front lenses and can be cleaned with a cloth or lens tissue wrapped around a finger and moistened with methanol. The 40X and 100X objectives require a little more care and examination with a magnifier is recommended.

To achieve the high degree of flatness obtained with the higher magnification objectives, it is necessary to use a small concave front lens of fairly short radius or curvature. The surface of this front lens can be readily cleaned with a toothpick covered with a cotton tip, or with a small cotton swab. Moisten the cotton with methanol and squeeze it almost dry. Wipe the front lens lightly without applying undue force or scrubbing action. Make sure that the cotton tip contacts the concave lens surface. Check the objective with a magnifier after cleaning.

Should you need to remove the viewing body, be careful not to accidentally touch the outer telelens surface on the underside of the body. Finger prints on this surface will reduce image clarity. This lens can be cleaned in the same manner as objectives and eyepieces.

### **6.2 Mechanical Maintenance**

Using the dust cover is the most important factor in keeping the microscope in good mechanical and physical condition.

The LEICA ATC™ 2000 has a durable finish. All non-optical parts of the instrument can be cleaned with methanol, naphtha or soap and water. Avoid the use of all other organic solvents when cleaning the instrument.

There is very little lubrication or routine maintenance required on the LEICA ATC™ 2000. The focusing mechanism does not require lubrication, nor does the mechanical stage. The only mechanism that requires periodic lubrication and cleaning is the substage condenser slide. This should be done by wiping away old lubricant with xylene and lubricating it with vaseline or a light grease.

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### 6.3 Lamp Replacement

**CAUTION:** The glass envelope of the lamp may be extremely hot. Do not attempt to change the lamp before it is completely cooled or without wearing adequate skin protection.

The only component in the LEICA ATC™ 2000 that may require periodic replacement is the illuminator lamp. The LEICA ATC™ 2000 utilizes a 6V, 9W halogen lamp (Part Number 631; Phillips Part Number 13529).

**CAUTION:** Lamp is not user serviceable. Please refer service to qualified service person.

- A. Before replacement, unplug the instrument.
- B. Unscrew the base assembly.
- C. After the lamp has cooled, loosen the spring retainer which holds the lamp in place and then carefully remove the lamp from its socket.
- D. Place the new lamp into the socket and reconnect the spring retainer over the lamp.
- E. Reattach the base assembly.
- F. Never operate the microscope illuminator unless the lamp window door is securely closed.

### 6.4 Electrical Considerations/Equipment Ratings

As a safety precaution, the power has been grounded to the microscope base. Never use a 2-prong adapter between the power cord and the power source; it will render the grounding feature ineffective.

**CAUTION:** Risk of electric shock – Removal of the bottom cover of the microscope exposes hazardous voltages, which if contacted, could cause injury or death. Please refer servicing to qualified personnel.

#### Supply Voltage and Range Supply Frequency

120V ~ ±10% - 60Hz

240V ~ ±10% - 50Hz

Input Power - 20W

#### Range of Environmental Conditions

- Altitude up to 2000m
- Indoor use only
- Temperatures 5°C to 40°C
- Maximum Relative Humidity - 80% for temperatures up to 31°C decreasing linearly to 50% relative humidity at 40°C
- Installation Category II (Overvoltage Category)
- Pollution degree 2

### 7.0 FUSE REPLACEMENT

**CAUTION:** Fuse Replacement - The fuse is not user serviceable. It must be replaced by qualified service personnel. The fuse can be replaced only with the same size, type and rating of fuse.

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*Printed in USA ©1996, Leica Inc., Buffalo NY USA 3/96 452RD-101 Rev. B*
*Leica Inc.*